

MD 100 Photometer

NEW

Precise Water
Analysis in High-
Quality Design

Small | Mobile | Rapid


▶ Please see pages 76 onwards for
tests, ranges and reagents



Highlights

- Waterproof*)
- Automatic Switch-Off
- Real-Time-Clock and Date
- Calibration Mode
- Backlit Display
- Storage Function
- Suitable for use as a testing instrument

*) as defined in IP 68, 1 hour at 0.1 meter

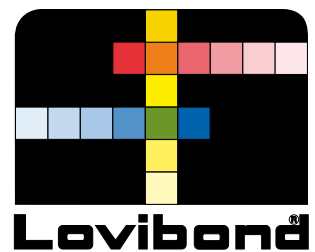

MD 100

Test	Code
Ammonia , tablet reagents 0.02 - 1.0 mg/l N	27 60 60
Ammonium , powder reagents 0.01 - 0.8 mg/l N	27 60 65
Ammonia, free powder reagents 0.01 - 0.5 mg/l N	27 60 70
Monochloramine 0.04 - 4.5 mg/l Cl ₂	
Chlorine , tablet reagents (OTZ) 0.01 - 6.0 mg/l Cl ₂ alternatively with liquid reagents 0.02 - 4 mg/l Cl ₂	27 60 00
Chlorine DUO , for 2 types of reagents 1) Tablet reagents 0.01 - 6.0 mg/l Cl ₂ 2) Powder reagents 0.02 - 2.0 mg/l Cl ₂ (ø 24 mm glass vial) 0.1 - 8.0 mg/l Cl ₂ (ø 10 mm multi vial-2)	27 60 20
Chlorine , powder reagents 0.02 - 2.0 mg/l Cl ₂ (ø 24 mm glass vial) 0.1 - 8.0 mg/l Cl ₂ (ø 10 mm multi vial-2)	27 60 30
Chlorine dioxide , tablet reagents 0.05 - 11 mg/l ClO ₂	27 60 35
Chlorine dioxide , powder reagents 0.05 - 3.8 mg/l ClO ₂	27 60 35
COD , tube tests 0 - 150 mg/l O ₂ (ø 16 mm) 0 - 1500 mg/l O ₂ (ø 16 mm) 0 - 15000 mg/l O ₂ (ø 16 mm)	27 61 20
Iron , tablet reagents 0.02 - 1.0 mg/l Fe	27 60 50

The MD 100 uses high quality interference filters with long-life LEDs as a light source without any moving parts in a transparency sample chamber.

The units supply accurate, reproducible results very quickly. Other major advantages include ease of operation, ergonomic design, compact dimensions and safe handling.

The calibration and software-based adjustment options mean that the MD 100 is also suitable for use as a testing instrument.

The tests are conducted using either Lovibond® tablet reagents with long-term stability and a guaranteed minimum 5 or 10 year shelf life, VARIO powder reagents or using liquid reagents.

Scroll Memory (SM)

To avoid unnecessary scrolling for the required test method, the instrument memorizes the last method used before switching off the instrument. When the instrument is switched on again, the scroll list comes up with the last used test method first.

MD 100

Test	Code
Iron TPTZ , powder reagents 0.02 - 1.8 mg/l Fe	27 60 55
Iron , powder reagents 0.02 - 3.0 mg/l Fe	27 60 56
Fluoride , without reagents 0.05 - 2.0 mg/l F ⁻	27 60 90
Copper , tablet reagents 0.05 - 5.0 mg/l Cu	27 60 80
Copper , powder reagents 0.05 - 5.0 mg/l Cu	27 60 85
Manganese LR , tablet reagents 0.2 - 4.0 mg/l Mn	27 61 00
Manganese LR , powder reagents 0.01 - 0.7 mg/l Mn	27 61 05
Manganese HR , powder reagents 0.1 - 18 mg/l Mn	27 61 06
Monochloramine powder reagents 0.04 - 4.5 mg/l Cl ₂	27 60 70
Phosphate , tablet reagents 0.05 - 4.0 mg/l PO ₄	27 60 40
Phosphate , powder reagents 0.06 - 2.5 mg/l PO ₄	27 60 45
Silica , tablet reagents 0.05 - 4.0 mg/l SiO ₂	27 61 10
Silica LR , powder reagents 0.1 - 1.6 mg/l SiO ₂	27 61 15
Silica HR , powder reagents 1 - 90 mg/l SiO ₂	27 61 16

Zero Setting (OTZ)

For certain versions of the instrument it is not necessary to zero the instrument each time. The zero setting is held in memory (One Time Zero - OTZ). The zero setting can be confirmed whenever it is useful.

➔ Please see pages 76 onwards for tests, ranges and reagents



MD 100 2in1

Test	Code
Chlorine, pH , tablet reagents (OTZ) 0.01 - 6.0 mg/l Cl ₂ ; 6.5 - 8.4 pH	27 80 20
Chlorine, pH , powder reagents for chlorine 0.02 - 2.0 mg/l Cl ₂ (ø 24 mm glass vial) 0.1 - 8.0 mg/l Cl ₂ (ø 10 mm multi vial-2) 6.5 - 8.4 pH	27 80 30

MD 100 3in1

Test	Code
Chlorine, pH, Cyanuric acid tablet reagents (OTZ) 0.01 - 6.0 mg/l Cl ₂ ; 6.5 - 8.4 pH 2 - 160 mg/l Cyanuric acid	27 80 10
Chlorine, pH, Alkalinity, total (M) tablet reagents (OTZ) 0.01 - 6.0 mg/l Cl ₂ ; 6.5 - 8.4 pH 5 - 200 mg/l CaCO ₃	27 80 60
Chlorine LR, Chlorine HR, Chlorine dioxide , tablet reagents 0.01 - 6.0 mg/l Cl ₂ 5 - 200 mg/l Cl ₂ (ø 16 mm round vial) 0.05 - 11 mg/l ClO ₂	27 80 00

MD 100 4in1

Test	Code
Chlorine, pH, Cyanuric acid, Alkalinity, total (M) tablet reagents (OTZ) 0.02 - 6.0 mg/l Cl ₂ ; 6.5 - 8.4 pH 2 - 160 mg/l Cyanursäure ; 5 - 200 mg/l CaCO ₃	27 80 70

MD 100 5in1

Test	Code
Chlorine, pH, Cyanuric acid, Alkalinity, total (M), Calcium hardness tablet reagents (OTZ) 0.02 - 6.0 mg/l Cl ₂ ; 6.5 - 8.4 pH 2 - 160 mg/l Cyanursäure ; 5 - 200 mg/l CaCO ₃ 0 - 500 mg/l CaCO ₃ (CaH)	27 80 80

MD 100 6in1

Test	Code
Chlorine, Bromine, pH, Cyanuric acid, Alkalinity, total (M), Calcium hardness tablet reagents (OTZ) 0.02 - 6.0 mg/l Cl ₂ ; 0.05 - 13 mg/l Br 6.5 - 8.4 pH ; 2 - 160 mg/l Cyanursäure 5 - 200 mg/l CaCO ₃ ; 0 - 500 mg/l CaCO ₃ (CaH)	27 80 90

Applications

- Water Treatment
- Waste/Drinking Water
- Cooling/Boiler Water
- Pool-Water
- Beverage Industry
- Laboratory and Field Testing

MD 100 Photometer



Calibration certificate

Beside the "Certificate of Compliance" which is supplied with the MD 100, calibration certificates are available at cost on request. Calibration certificates are individually supplied per instrument and per method.

The Calibration Certificate has to be ordered together with the new instrument and cannot be delivered at a later stage.

N.I.S.T Traceability

The instrument has a factory calibration, which is related to international standards which are not N.I.S.T traceable. The instrument may be calibrated by the user in a "user calibration mode" with N.I.S.T traceable standards.

(N.I.S.T. = National Institute of Standards and Technology)

Delivery Content

- Instrument in carrying case
- 4 micro batteries (AAA)
- 3 round vials (glass) with lids
- Tablet reagents and/or liquid reagents or VARIO Powder reagent
- Guarantee sheet
- Certificate (COC)
- Instruction Manual

You can find updated information on parameters and measuring ranges on our website at www.tintometer.com

➔ Please see pages 76 onwards for tests, ranges and reagents

Technical Data

Optics LEDs, interference filters (IF) and photo sensor in transparent sample chamber. Depending on the version, up to 3 different interference filters are used. Wavelength specifications of interference filters:
 430 nm $\Delta\lambda = 5$ nm
 530 nm $\Delta\lambda = 5$ nm
 560 nm $\Delta\lambda = 5$ nm
 580 nm $\Delta\lambda = 5$ nm
 610 nm $\Delta\lambda = 6$ nm
 660 nm $\Delta\lambda = 5$ nm

Wavelength Accuracy	± 1 nm
Photometric Accuracy⁴⁾	3% FS (T = 20°C – 25°C)
Photometric Resolution	0.01 A
Power Supply	4 micro batteries (AAA), capacity approx. 17 hours or 5000 tests
Auto - OFF	automatic switch-off
Display	backlit LCD (on keypress)
Storage	internal ring memory for 16 data sets
Interfaces	Infrared interface for test data transfer
Additional feature	real time clock and date
Calibration	factory calibration and user calibration. Reset to factory calibration possible
Dimensions	155 x 75 x 35 mm (L x W x H)
Weight	basic unit approx. 260 g
Environmental conditions	Temperature: 5 – 40 °C rel. humidity: 30 – 90% (non condensing)

CE-Conformity

⁴⁾ tested with standard solutions

Accessories

Item	Code
Set of 12 round vials with cap Height 48 mm, Ø 24 mm	19 76 20
Set of 10 round vials with cap Height 90 mm, Ø 16 mm	19 76 65
Adapter for round vials Ø 16 mm	19 80 22 20
Set of 12 plastic vials (PC), with lid "Multi"-Type 2	19 76 00
Cleaning cloth for vials	19 76 35
Measuring beaker, volume 100 ml	38 48 01
Plastic funnel with handle	47 10 07
Plastic stirring rod, 13 cm length	36 41 00
Plastic stirring rod, 10 cm length	36 41 09
4 micro batteries (AAA)	19 50 026
Infra-red data transfer modul IRiM	21 40 50



Data transfer

The optional available IRiM (infra-red interface modul) uses modern infra-red technology to transmit measurement data from the MD 100 photometer to one of 3 optional interfaces. These interfaces can be used to connect to a PC, a USB printer¹⁾ or alternative a serial printer²⁾.

The unit is supplied complete with data logging software providing easy and rapid transfer of data to the PC. As an option the data can be saved as an Excel sheet or a .txt file.

Measurement data can quickly be printed out, using a specified¹⁾ USB or alternative a printer with a serial plug-in connected to the IRiM.

¹⁾ USB printer: HP Deskjet 6940 ; ²⁾ each ASCII printer

Verification Standard Kit

The Verification standard kit for the MD 100 is designed to reassure the user about the accuracy and the reliability of the results.

The shelf life of the Verification standard kit is two years from the date of production, provided that storage and use are in accordance with the instructions provided.

Verification Standard Kit 21 56 70



Reference Standard Kit

The reference standards are designed to check the accuracy and the reliability of the results.

The shelf life of reference standards is two years from the date of production, provided that storage and use are in accordance with the instructions provided.

Kit Chlorine for instruments with tablet / liquid reagent 0.2* and 1.0* mg/l 27 56 50

Kit Chlorine for instruments with tablet / liquid reagent 0.5* and 2.0* mg/l 27 56 55

Kit Chlorine for instruments with tablet / liquid reagent 1.0* and 4.0* mg/l 27 56 56

Kit Chlorine for instruments with powder reagent (VARIO) 0.2* and 1.0* mg/l 27 56 60

Kit pH for instruments with tablet / liquid reagent 7,45* pH 27 56 70

* Approximate figure, actual figure specified in Certificate of Analysis